

Near Infrared Spectroscopy: The Future Waves

Edited by
A.M.C. Davies and
Phil Williams

The Proceedings of the 7th International Conference on
Near Infrared Spectroscopy, Montréal, Canada,
6-11 August 1995

NIR *Publications*

Published by:

NIR Publications, 6 Charlton Mill, Charlton, Chichester, West Sussex PO18 0HY, United Kingdom, Tel: +44(0)1243-811334, Fax: +44(0)1243-811711.

ISBN: 0 9528666 0 9

British Library Cataloguing-in-Publication Data

A catalogue record for this book is available from the British Library


© NIR Publications 1996

All rights reserved. Apart from any fair dealing for the purposes of research or private study, or criticism or review, as permitted under the UK's Copyright, Designs and Patents Act 1988, this publication may be reproduced, stored or transmitted, in any form or by any means, only with the prior permission in writing of the publishers, or in the case of reprographic reproduction, in accordance with the terms of licences issued by the Copyright Licensing Agency. Enquiries concerning reproduction outside those terms should be sent to the publishers at the address above.


Typeset by IM Publications, Chichester, UK.

Printed and bound in the UK by H.W. Charlesworth Ltd, Huddersfield.

Ne

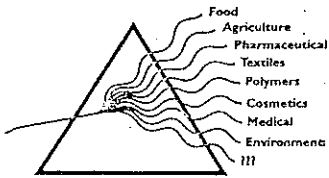

Near Infrared Spectroscopy: The Future Waves

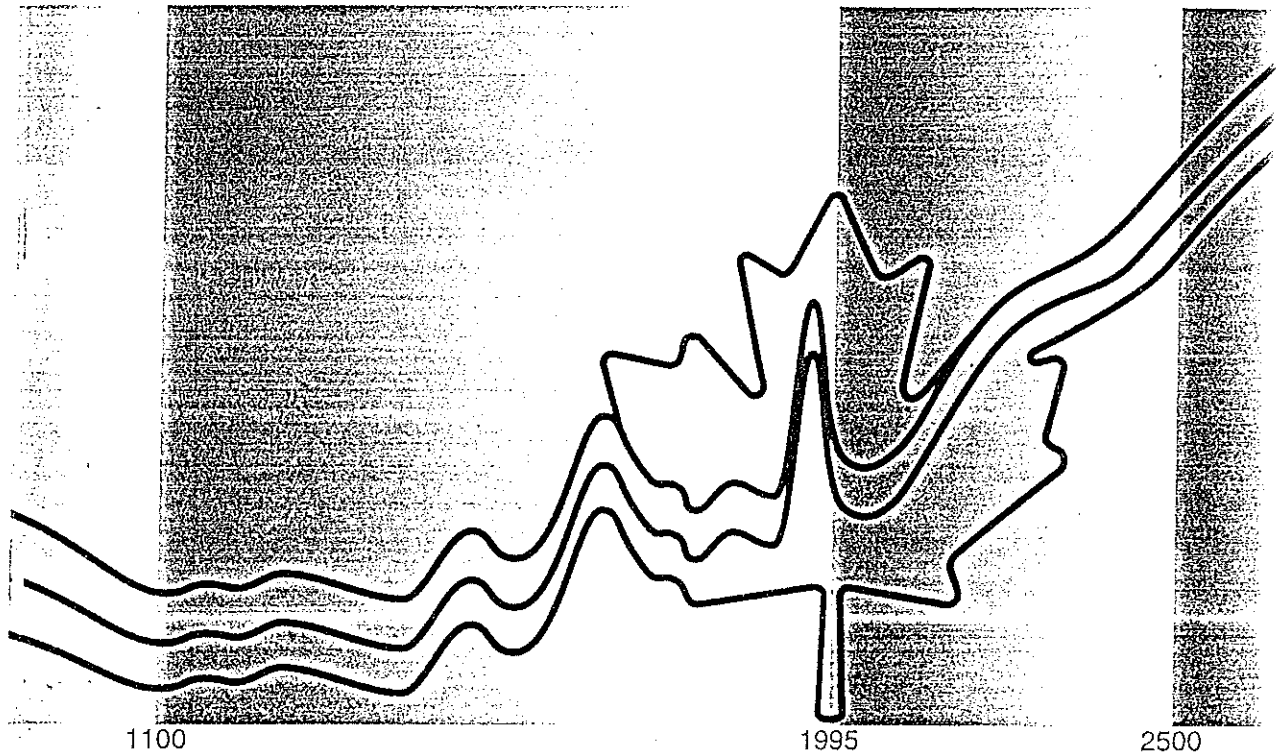

Edited by
A.M.C. Davies and Phil Williams

NIR Publications

Study on yeast enzymatic activities by near infrared spectroscopy and radiorespirometry

S.H. Yan, M. Maudoux, S. Collin, J-L. Van Haecht and J. Mayaudon

Catholic University of Louvain, Brewing and Food Industries, Place Croix du Sud 2-7, INBR Unit, 1348 Louvain-la-Neuve, Belgium.

Introduction

In the field of yeast fermentation, optimisation of biotechnological processes requires the application of a rapid and accurate method which offers the possibility of a feedback quality control.

The reliability of NIR spectroscopy was supervised by using a radiorespirometer fitted with an ionisation chamber.¹ Our aim is to study the applicability of near infrared (NIR) spectroscopy to measure the yeast enzymatic activities during the fermentation.

Material and methods

Yeast strain: *Saccharomyces cerevisiae* BRAS 268 was used in top fermentation (10°P wort). Storage was at 4°C on wort-agar slants. During fermentation, the yield of yeast cells was determined by dry weight expressed in mg 100 mL⁻¹ obtained at 50°C, *in vacuo*. The carbohydrate content was determined by the EBC method and expressed in g 100g⁻¹.

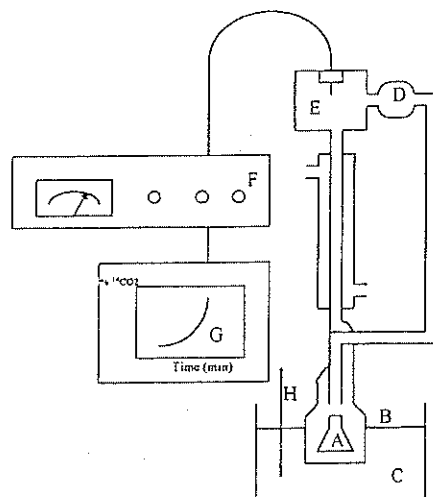


Figure 1. Scheme of radiorespirometry.¹

The production of yeast acetone powders was described elsewhere.² NIR spectra were recorded every 2 nm from 1100 to 2500 nm by using a spectrometer, PS Co 6250, working in a single beam mode. The active water content was determined by NIR spectroscopy and calculated by absorbance difference ($\log 1/R$)—A: between 1938 and 1858 nm and —B: between 1938 and 2020 nm.

The enzymatic activities were obtained by radiorespirometry (Figure 1). The data were expressed in % of maximum degradation of glucose-U-¹⁴C and initial velocity (V) of glucose-U-¹⁴C degradation min^{-1} .

Results and discussion

The transformation of glucose by yeasts during fermentation has been carried out under proliferating conditions. Measurements were made at different time plots varying from 5 to 200 hour incubation. Fermentation profiles of the various carbohydrate contents and yeast production are presented in Figure 2.

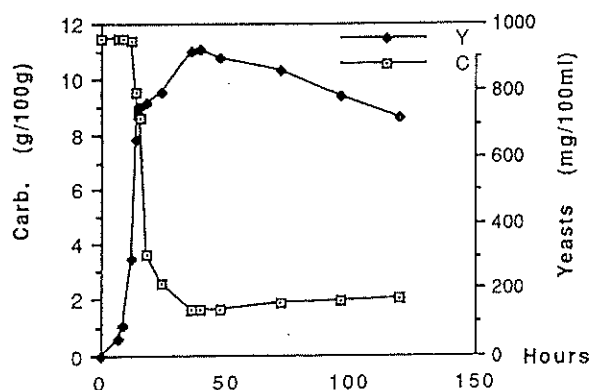


Figure 2. Carbohydrate content and yeast production.

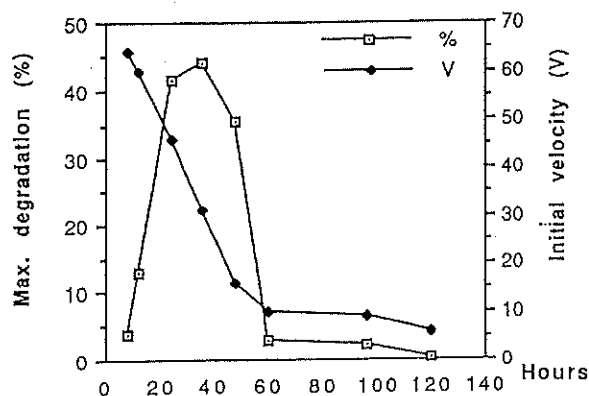


Figure 3. Initial velocity and degradation of glucose C14.

ra were recorded
in a single beam
ulated by absor-
38 and 2020 nm.
The data were
V) of glucose-U-

The results showed a linear relationship between, on the one hand, the breakdown of glucose-U-¹⁴C and the enzymatic activities (Figure 3) and on the other hand, active water content in living yeast cells (Figure 4 and Figure 5).

carried out under
g from 5 to 200
yeast production

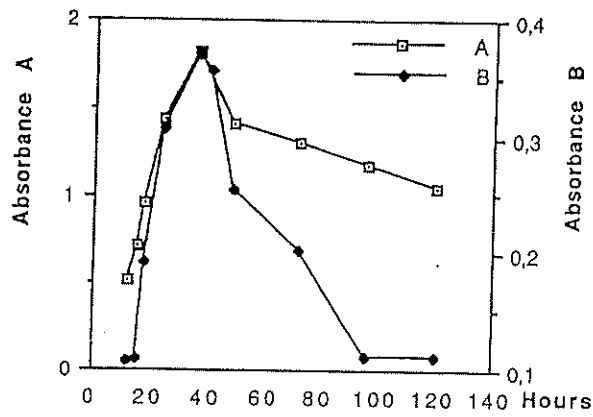


Figure 4. Active water contents during yeast fermentation.

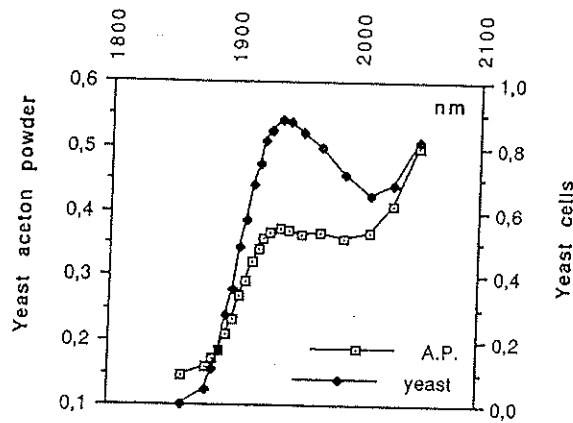


Figure 5. NIR spectra of yeast and yeast acetone powder.

Table 1: Maximum degradation (%¹⁴CO₂) in different reaction systems.

Systems with	glucose 2- ¹⁴ C	glucose 3,4- ¹⁴ C
ATP = 0.01 mg	25	40
ATP = 0	0	0

Glucose 2-¹⁴C or glucose 3,4-¹⁴C, 1 micro Ci/nM; yeast acetone powder 50 mg; 1/15 M phosphate buffer 1 mL, pH 4.6; Temp. 28°C.

Let us note, by using glucose-U- ^{14}C and glucose-2- ^{14}C (instead of glucose-U- ^{14}C) that with yeast acetone powders, enzymes of Embden-Meyerhof-Parnas and tricarboxylic acid pathways were active only if ATP traces were added in the mixture reaction systems as described in Table 1.

Conclusion

These results show that NIR spectroscopy, with the help of radiorespirometry, could provide the possibility to rapidly follow any yeast fermentation process for qualitative and quantitative purposes.

References

1. J. Mayaudon, *Soil Biochemistry* **2**, 202 (1971).
2. I.C. Gunsalus, *Methods in Enzymology* **1**, 51 (1955).

Applica spectro analysis

I.K. Brookes,

Pfizer International

Introduction

Since 1988 F
fermentation pro.
International NIF
recovery and the

Experimenta

Safety

NIR analysis
Appropriate prec

Apparatus

Laboratory i
mator (NIRSystem
module for all at
reflectance fibre
NIRSystems "Ir
trophotometer C
with the instrum

Procedure for at-

Calibration
methods and co
be representati
formulation and
each parameter
calibration for :